

**GERSTEL**

# AppNote 4/2010

## Automated QuEChERS Extraction for the Confirmation of Pesticide Residues in Foods using LC/MS/MS

Virgil Settle, Fred Foster  
*Gerstel, Inc., 701 Digital Dr. Suite J,  
Linthicum, MD 21090, USA*

Paul Roberts  
*Anatune, Ltd., Broadway House, 148-152 St. Neots Road, Hardwick,  
Cambridge, CB23 7QJ, UK*

Peter Stone  
*Agilent Technologies, 5301 Stevens Creek Blvd., Santa Clara,  
CA 95051-7201, USA*

Joan Stevens  
*Agilent Technologies, 2850 Centerville Rd.,  
Wilmington, DE 19808, USA*

Jon Wong, Kai Zhang  
*US FDA, 5100 Paint Branch Parkway, HFS-3 36, College Park,  
MD 20740, USA*

### **KEYWORDS**

QuEChERS, LC/MS/MS, Sample Preparation, Lab Automation

### **ABSTRACT**

QuEChERS (quick, easy, cheap, effective, rugged, and safe) sample preparation methods were developed to help monitor pesticides in a range of food samples [1]. These methods are quite labor intensive, however, since they include several manual steps, such as shaking, centrifugation, and dispersive SPE. If dispersive SPE clean up of QuEChERS type extracts could be automated, laboratory productivity for monitoring pesticide residues in food could be improved significantly. In the work presented here, an automated dispersive SPE clean-up method for QuEChERS extracts was developed and combined with LC/MS/MS analysis of the cleaned

extracts. Extracts were prepared with commercially available kits from Agilent Technologies. Automation was achieved using a GERSTEL MPS XL PrepStation configured with an Anatune CF-100 centrifuge. Analytical methodology for confirming the presence of a variety of pesticides in a range food samples was developed using an Agilent G6460A Triple Quadrupole Mass Spectrometer. The sensitivity and selectivity of LC/MS/MS enable sufficiently low limits of determination to meet acceptance criteria for reporting the maximum residue levels (MRLs) that are established by regulatory agencies.

The ability to automate the dispersive SPE clean-up of QuEChERS extracts combined with introduction of the cleaned extract directly to the LC/MS/MS system results in improved laboratory productivity by streamlining the complete analytical process.

## INTRODUCTION

Automated QuEChERS clean-up of fruit and vegetable extracts combined with LC/MS/MS determination of pesticides was previously reported [2]. Recovery and reproducibility data based upon matrices spiked at concentrations of 10 ng/g showed that the system used was suitable for automating the QuEChERS dispersive SPE clean-up.

Recently, interest was expressed in an automated method to determine the presence of pesticide residues in botanical samples. Some botanical matrices have physical properties that pose a real challenge. When using the automated QuEChERS clean-up procedure, such matrices can make it difficult to reach the low limits of determination required in order to meet acceptance criteria for reporting the maximum residue levels (MRLs) as established by regulatory agencies.

This study focuses on the automated extraction used in the second step of the QuEChERS procedure, which is followed by LC/MS/MS analysis of the extract. The aim is to provide high throughput analysis for the confirmation of pesticide residues in botanical matrices. Automated QuEChERS extractions were performed by the GERSTEL MPS XL autosampler using Agilent's SampliQ QuEChERS dispersive SPE sorbent blend suggested for fatty matrices. The resulting extracts from the automated QuEChERS process were introduced into an Agilent 6460 LC/MS/MS instrument.

## EXPERIMENTAL

*Materials.* Stock solutions containing the pesticide compounds listed in Table 1 in acetonitrile were prepared and provided by the FDA. Calibration standards and matrix matched standards were prepared by making appropriate dilutions of the pesticide stock solutions using mobile phase, blank hop extract, or blank ginseng extract resulting in the following concentrations: 0.5, 1, 2, 5, 10, 20, 50, 100, 200, 500, and 1000 ng/mL.

Crude acetonitrile extracts of pesticide-fortified samples, incurred samples, and blank matrix samples based on both Hops and Ginseng root were prepared and provided by the FDA. These samples were generated using QuEChERS extraction salts for the EN Method (Agilent p/n: 5982-5650) and the recommended sample preparation method supplied with the salts [3].

*Instrumentation.* All automated PrepSequences were performed using an MPS 2XL dual head, multi-purpose sampler configured for automated QuEChERS-LC/MS/MS analysis as shown in Figure 1.



**Figure 1.** MPS 2XL dual head, multi-purpose sampler configured for automated QuEChERS-LC/MS/MS analysis.

**Table 1.** 200+ pesticides monitored using the automated QuEChERS-LC/MS/MS method.

|                     |                     |                    |                      |                    |
|---------------------|---------------------|--------------------|----------------------|--------------------|
| 3-Hydroxycarbofuran | Acephate            | Acetamiprid        | Acibenzolar-S-methyl | Alanycarb          |
| Aldicarb            | Aldicarb sulfone    | Aldicarb sulfoxid  | Aspon                | Avermectin B1a     |
| Avermectin B1b      | Azadirachtin        | Azoxystrobin       | Benalaxyl            | Bendiocarb         |
| Benfuracarb         | Benoxacor           | Benthiavalicarb    | Benzoximate          | Bifenazate         |
| Bifenthrin          | Bitertanol          | Boscalid           | Bromuconazole-1      | Bromuconazole-2    |
| Bupirimate          | Buprofezin          | Butafenacil        | Butocarboxym         | Butoxycarboxim     |
| Cadusafos           | Carbaryl            | Carbendazim        | Carbetamid           | Carbofuran         |
| Carboxine           | Carfentrazone-ethyl | Chlordimeform      | Chlorfenvinphos-beta | Chlorfluazuron     |
| Chlorotoluron       | Chloroxuron         | Clethodim          | Clofentezine         | Clothianidin       |
| Coumaphos           | Cumyluron           | Cyanazine          | Cyanophos            | Cyazofamid         |
| Cycluron            | Cymoxanil           | Cyproconazole      | Cyprodinil           | Cyromazine         |
| d10-Diazinon        | d6-Dichlorvos       | d6-Dimethoate      | d6-Diuron            | d6-Linuron         |
| d6-Malathion        | Daimuron            | Dazomet            | Deltamethrin         | Diazinon           |
| Dichlorvos          | Dicrotophos         | Diethofencarb      | Difenoconazol        | Diflubenzuron      |
| Dimethenamid        | Dimethoat           | Dimethomorph A     | Dimethomorph B       | Dimoxystrobin      |
| Diniconazole        | Dinotefuran         | Dioxacarb          | Disulfoton           | Dithiopyr          |
| Diuron              | Dodemorph 1         | Dodemorph 2        | E-Fenpyroximate      | Emamectin B1a      |
| Emamectin B1b       | Epoconazole         | Eprinomectin B1a   | EPTC                 | Esprocarb          |
| Ethidimuron         | Ethiofencarb        | Ethion             | Ethiprole            | Ethirimol          |
| Ethofumesate        | Ethoprop            | Etobenzanid        | Etofenprox           | Etoazole           |
| Famoxadone          | Fenamidone          | Fenarimol          | Fenazaquin           | Fenbuconazol       |
| Fenhexamid          | Fenoxanil           | Fenoxycarb         | Fenpropathrin        | Fenpropimorph      |
| Fenuron             | Flonicamid          | Flucarbazone       | Fludioxinil          | Flufenacet         |
| Flufenoxuron        | Flumetsulam         | Flumioxazin        | Fluometuron          | Fluquinconazole    |
| Flusilazol          | Fluthiacet-methyl   | Flutolanil         | Flutriafol           | Forchlorfenuron    |
| Formetanate         | Fuberidazole        | Furalaxyl          | Furathiocarb         | Heptenophos        |
| Hexaconazol         | Hexaflumuron        | Hexythiazox        | Hydramethylnon       | Imazail            |
| Imazapyr            | Imibenconazole      | Imidacloprid       | Indanofan            | Indoxacarb         |
| Ipconazole          | Iprovalicarb        | Isocarbamid        | Isofenfos            | Isopropalin        |
| Isoproturon         | Isoxaben            | Isoxaflutole       | Kresoxim-methyl      | Lactofen           |
| Leptophos           | Linuron             | Lufenuron          | Mandipropamid        | Mefenazet          |
| Mepanipyrim         | Mepronil            | Metalaxyl          | Metconazole          | Methabenzthiazuron |
| Methamidophos       | Methiocarb          | Methomyl           | Methoprotryne        | Methoxifenozyd     |
| Metobromuron        | Metribuzin          | Mevinphos          | Mexacarbate          | Molinate           |
| Monocrotophos       | Monolinuron         | Moxidectin         | Myclobutanil         | Neburon            |
| Nitenpyram          | Norflurazon         | Novaluron          | Nuarimol             | Omethoate          |
| Oxadixyl            | Oxamyl              | Paclobutrazol      | Penconazole          | Pencycuron         |
| Phenmedipham        | Picoxystrobin       | Piperonyl butoxide | Pirimicarb           | Prochloraz         |
| Promecarb           | Prometon            | Prometryn          | Propachlor           | Propamocarb        |
| Propargite          | Propazine           | Propham            | Propiconazole        | Propoxur           |
| Pymetrozine         | Pyracarbolid        | Pyraclostrobin     | Pyridaben            | Pyrimethanil       |
| Pyriproxyfen        | Quinoxifen          | Rotenone           | Sebuthylazine        | Secbumeton         |
| Siduron             | Simazine            | Simetryn           | Spinosyn A           | Spinosyn D         |
| Spirodiclofen       | Spiromesifen        | Spiroxamin         | Sulfentrazone        | Tebuconazole       |
| Tebufenozide        | Tebufenpyrad        | Tebuthiuron        | Teflubenzuron        | Temephos           |
| Terbumeton          | Terbutryn           | Terbutylazine      | Tetraconazole        | Tetramethrin cis   |
| Thiabendazole       | Thiacloprid         | Thiametoxam        | Thiazopyr            | Thidiazuron        |
| Thiobencarb         | Thiofanox           | Thiophanate-methyl | Triadimefon          | Triadimenol        |
| Trichlamide         | Trichlorfon         | Tricyclazole       | Trifloxystrobin      | Triflumizole       |
| Triflumuron         | Triticonazole       | Uniconazole        | Vamidotion           | Zoxamide           |

### QuEChERS acetonitrile extract pretreatment.

- Pipette 1mL of the acetonitrile extract obtained following the 1st centrifugation step of the QuEChERS sample preparation method into a 2 mL glass autosampler vial containing the sorbent blend from Agilent's SampliQ QuEChERS dispersive SPE kit for fatty samples, AOAC (p/n: 5982-5122).
- Place the sample onto a tray on the GERSTEL MPS 2XL dual head, multi-purpose sampler configured for automated QuEChERS-LC/MS/MS.

The automated QuEChERS extraction consisted of the following steps:

### Automated QuEChERS Prep Sequence.

- Agitate the sample vial for 1 minute using the Anatune CF-100 centrifuge.
- Centrifuge the sample vial at 575 g for 3 minutes using the Anatune CF-100 centrifuge.
- Filter 500 µL of the resulting supernatant through a 0.45 µm GERSTEL format syringe filter.
- Combine 100 µL of the resulting filtrate with 400 µL of mobile phase A in a clean 2 mL vial.
- Agitate the sample vial using the Anatune CF-100 centrifuge for 30 seconds.
- Inject 2 µL into the LC/MS/MS system.

Figure 2 shows a Prep Sequence used to perform an automated QuEChERS-LC/MS/MS run.

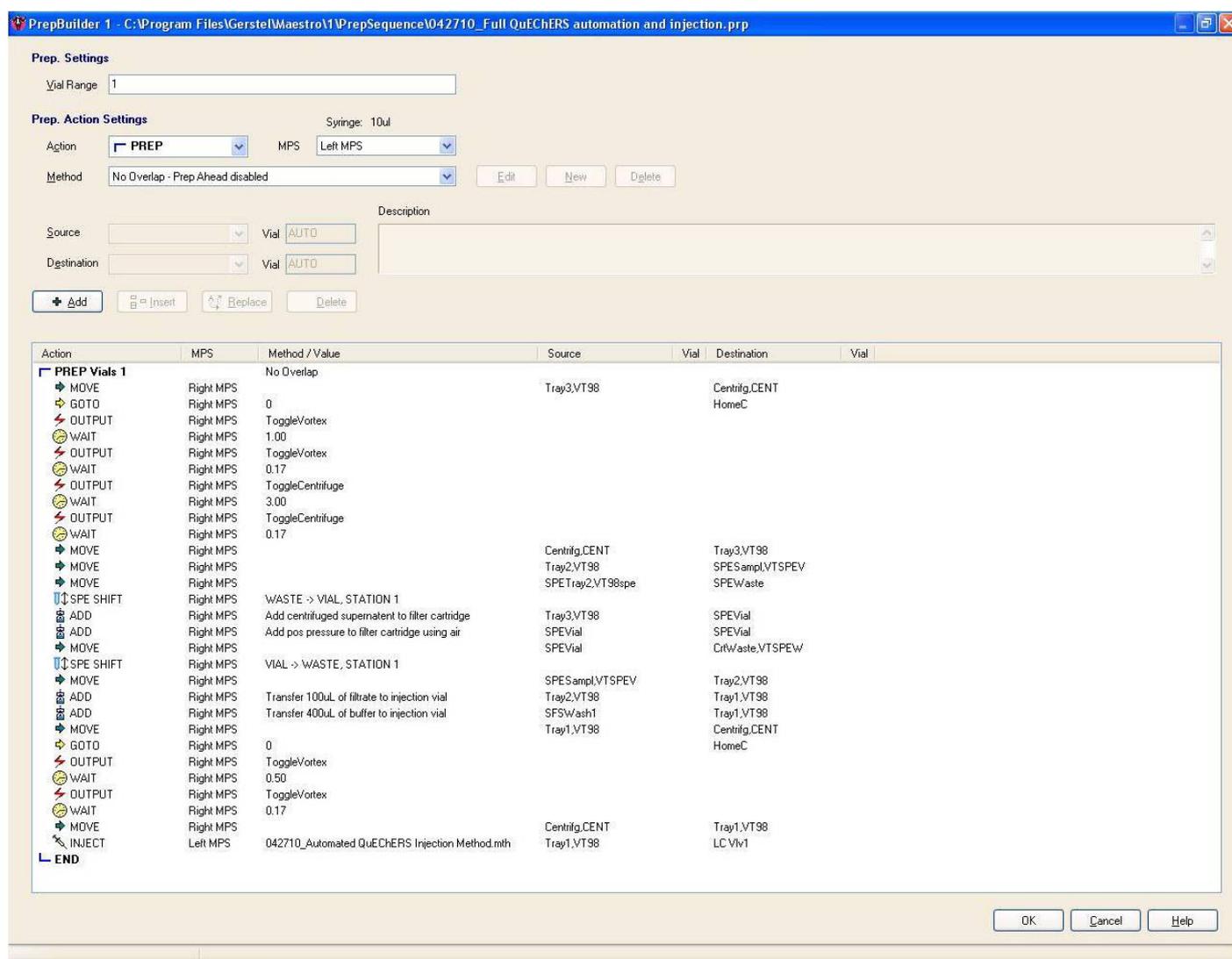


Figure 2. Example Prep Sequence for Automated QuEChERS LCMSMS Analysis.

Preparation of all standards was automated using the MPS 2XL dual head, multipurpose sampler configured for automated QuEChERS-LC/MS/MS analysis as follows:

*Preparation of Calibration Standards and Matrix Matched Standards.*

- Transfer 100  $\mu$ L of previously extracted matrix blank or 100 % acetonitrile into an empty 2 mL autosampler vial.
- Transfer 250  $\mu$ L of mobile phase A into the vial.
- Transfer 150  $\mu$ L of the respective standard stock solution into the vial.
- Agitate the vial using the Anatune CF-100 centrifuge for 30 seconds.

All analyses were performed using an Agilent 1290 HPLC, an Agilent 6460 Triple Quadrupole Mass Spectrometer with electrospray source and Jet Stream Option and a GERSTEL MPS 2XL autosampler configured with Active WashStation. Sample injections were made using a 6 port (0.25 mm) Cheminert C2V injection valve fitted with a 2  $\mu$ L stainless steel sample loop.

*Analysis conditions LC.*

Mobile Phase: A - 5 mM ammonium formate  
in water with 0.01 % formic acid  
B - 0.01 % formic acid  
in acetonitrile

Gradient: Initial 94 % A / 6 % B  
0.3 min 94 % A / 6 % B  
14 min 5 % A / 95 % B  
17 min 5 % A / 95 % B

Pressure: 600 bar  
Flowrate: 500  $\mu$ L/min  
Runtime: 17 min  
Post time: 2.5 min  
Column: 2.1 mm x 100 mm, 1.8  $\mu$ m,  
Zorbax Eclipse Plus C18 RRHT  
(Agilent)  
Oven: 55°C  
Inj. volume: 2  $\mu$ L


*Analysis conditions MS.*

Operation: Electrospray positive mode  
(Jet Stream)  
Time Filter Width: 0.04 min  
Scan Type: Dynamic MRM  
Delta EMV: 0 V  
Cycle Time: 660 ms  
Gas Temperature: 225 °C  
Gas Flow (N2): 10 L/min  
Nebulizer pressure: 25 psi  
Sheath Gas (N2): 350 °C  
11 L/min  
Capillary voltage: 4500 V  
Nozzle Voltage: 500 V

The mass spectrometer acquisition parameters and respective quantifier/qualifier ion transitions were chosen using the pesticide database option available for the MassHunter B.03.01 software. Table 1 provides a list of the more than 200 pesticides that were monitored using this single LC/MS/MS method. A retention time window value of 0.5 minute was used for each positive ion transition being monitored during the course of the dynamic MRM experiment.


**RESULTS AND DISCUSSION**

Figure 3 shows a representative overlay mass chromatogram resulting from a QuEChERS extract of a pesticide-fortified hop sample. More than 200 different pesticides were successfully determined in this botanical matrix using the automated QuEChERS-LC/MS/MS method.

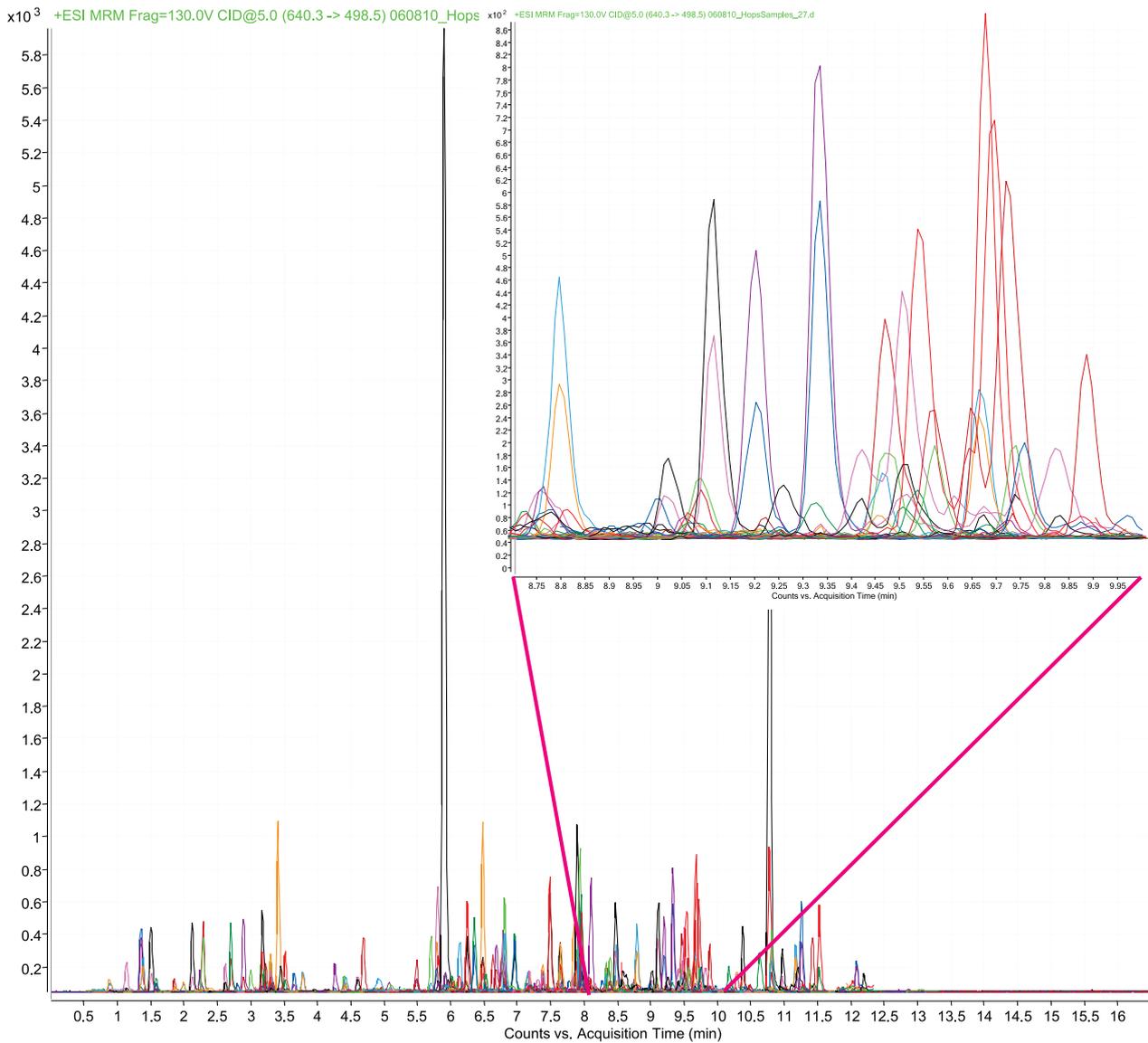


**Figure 3.** Representative Mass Chromatograms for low QC sample.

Figure 4 shows the a representative calibration curve resulting from automated preparation of neat standards. The calibration curves were shown to be linear from at least 1.00 to 200 ppb for the pesticides monitored, using a linear, 1/x regression method.

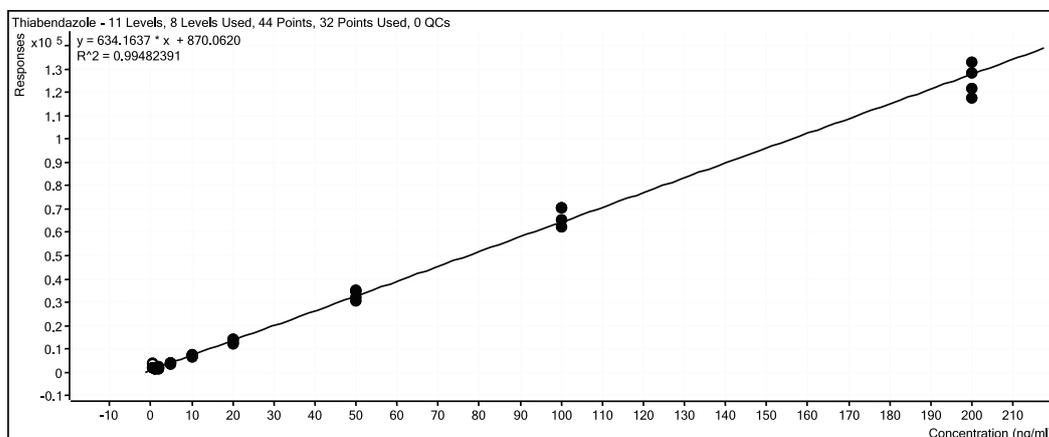


**Figure 4.** Representative calibration curve from automated neat standard preparation: thiabendazole.


Figures 5 through 7 show representative overlay mass chromatograms of neat, hop matrix matched, and ginseng matrix matched calibration standards respectively, all at a concentration of 10 ppb. The calibration standards were prepared automatically by the MPS 2XL.

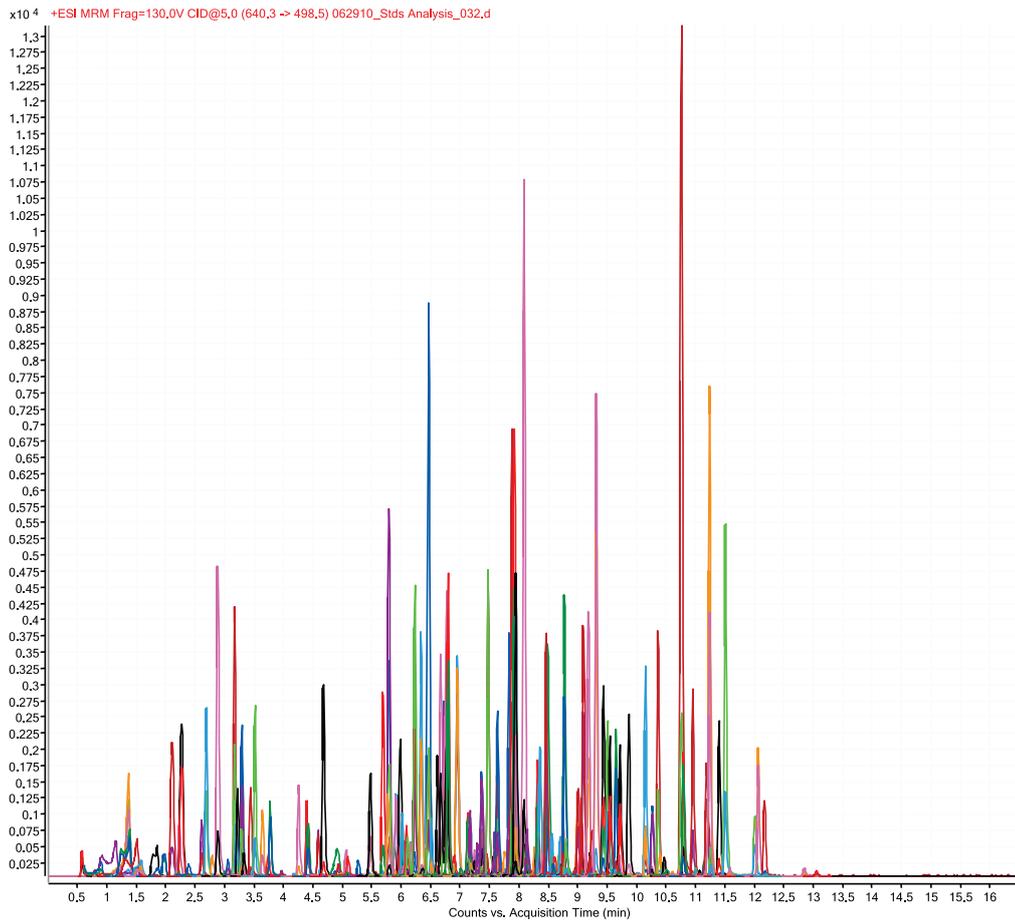

**Figure 5.** Representative overlay mass chromatogram for a 10 ppb neat standard.

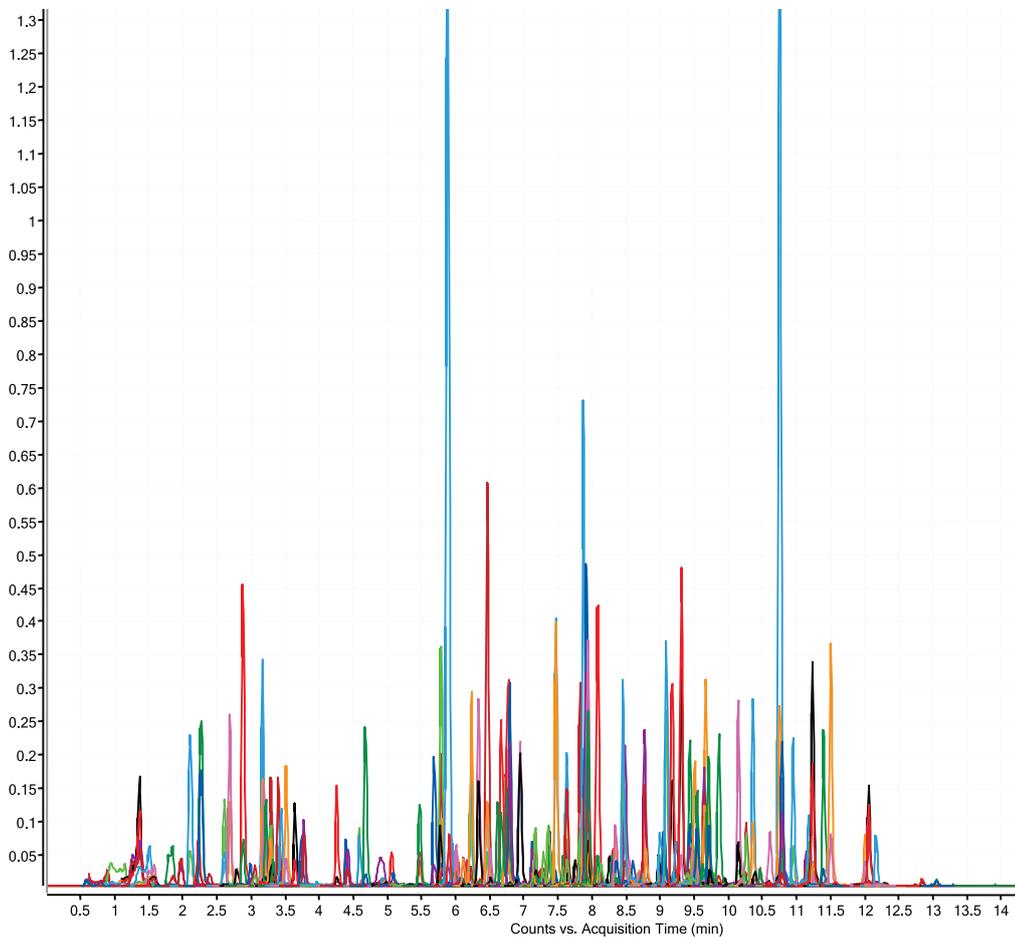

**Figure 6.** Representative overlay mass chromatogram for a 10 ppb hops matrix matched standard.

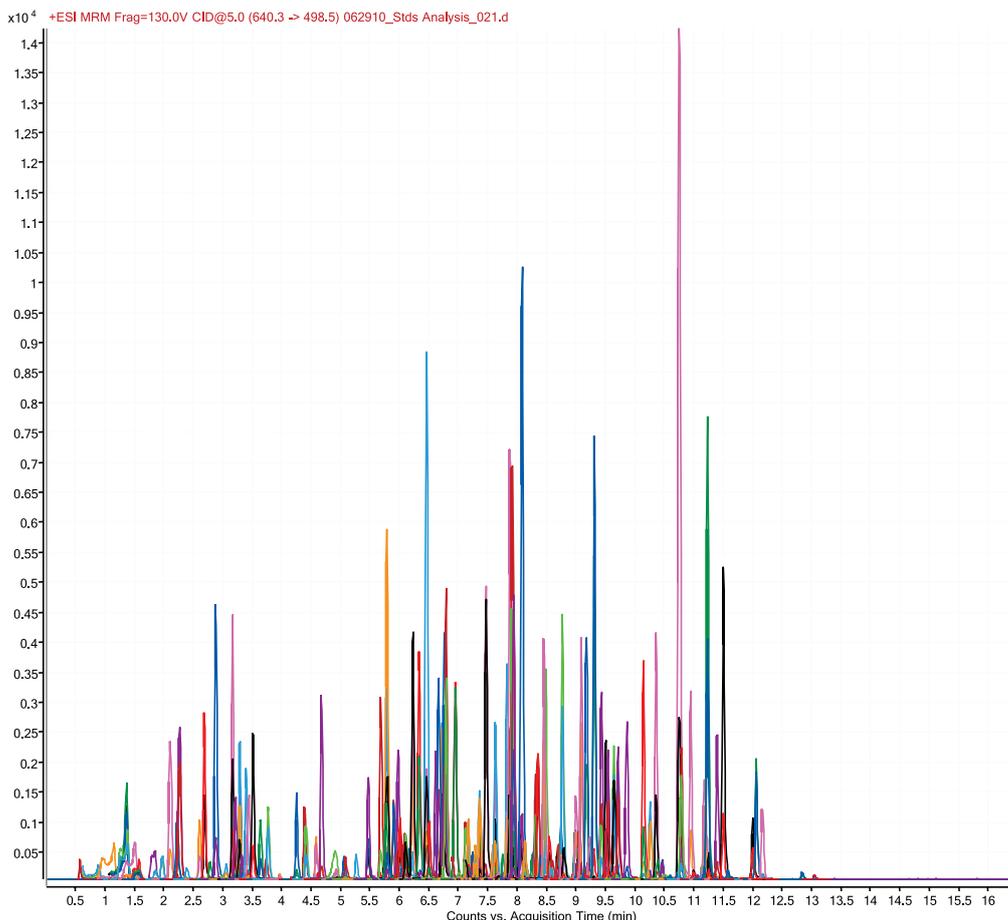

**Figure 7.** Representative overlay mass chromatogram for a 10 ppb ginseng matrix matched standard.

The total time required per sample to perform the QuEChERS extraction was 15 minutes. This was shorter than the LC-MS/MS analysis run, enabling the MPS system to complete preparation of the next sample during the LC-MS/MS run. Sample preparation and analysis are easily synchronized using the MAESTRO software providing “just in time” sample preparation and maximizing sample throughput.

## CONCLUSIONS

In this study, we were able to demonstrate:

- Successful monitoring of more than 200 pesticides in botanical matrix samples using an automated QuEChERS extraction method coupled to LC/MS/MS analysis using the Agilent 6460 Triple Quadrupole Mass Spectrometer.
- Automation of both the QuEChERS method and the preparation of standards using the GERSTEL MPS 2XL dual head robotic sampler.
- The “just-in-time” sample preparation capability included in the MAESTRO software enables highly efficient QuEChERS extraction and analysis.

## REFERENCES

- [1] Fully Automated QuEChERS clean-up and LC/MS-QQQ analysis of pesticides in Fruits and Vegetables. Anatune Chromatography Technical Note, No. AS90, Paul H. Roberts, Anatune Ltd., 2009.
- [2] Determination of pesticide residues in foods by acetonitrile extraction and partitioning with magnesium sulfate: collaborative study. Lehotay, S.J., J. AOAC Int. 90, 485 (2007)
- [3] Agilent SampliQ QuEChERS Kit brochure. Document number 5990-3562EN [www.chem.agilent.com](http://www.chem.agilent.com)





### GERSTEL GmbH & Co. KG

Eberhard-Gerstel-Platz 1  
45473 Mülheim an der Ruhr  
Germany

-  +49 (0) 208 - 7 65 03-0
-  +49 (0) 208 - 7 65 03 33
-  gerstel@gerstel.com
-  www.gerstel.com

## GERSTEL Worldwide

### GERSTEL, Inc.

701 Digital Drive, Suite J  
Linthicum, MD 21090  
USA

-  +1 (410) 247 5885
-  +1 (410) 247 5887
-  sales@gerstelus.com
-  www.gerstelus.com

### GERSTEL AG

Wassergrabe 27  
CH-6210 Sursee  
Switzerland

-  +41 (41) 9 21 97 23
-  gerstelag@ch.gerstel.com
-  www.gerstel.ch

### GERSTEL K.K.

1-3-1 Nakane, Meguro-ku  
Tokyo 152-0031  
SMBC Toritsudai Ekimae Bldg 4F  
Japan

-  +81 3 5731 5321
-  +81 3 5731 5322
-  info@gerstel.co.jp
-  www.gerstel.co.jp

### GERSTEL LLP

Level 25, North Tower  
One Raffles Quay  
Singapore 048583

-  +65 6622 5486
-  +65 6622 5999
-  SEA@gerstel.com
-  www.gerstel.com

### GERSTEL Brasil

Av. Pascoal da Rocha Falcão, 367  
04785-000 São Paulo - SP Brasil

-  +55 (11)5665-8931
-  +55 (11)5666-9084
-  gerstel-brasil@gerstel.com
-  www.gerstel.com.br



Awarded for the  
active pursuit of  
environmental sustainability

Information, descriptions and specifications in this Publication are subject to change without notice. GERSTEL, GRAPHPACK and TWISTER are registered trademarks of GERSTEL GmbH & Co. KG.

© Copyright by GERSTEL GmbH & Co. KG